

Athena

Multifunctional chair



SEATING

Frequent Usage

Chemotherapy
Infusions Therapy
Day Hospital
Blood donor
Blood samples
Emergency
Dialysis
Patient transport
Dermatology
Cosmetics
Multifunctional



WISDOM
that creates WELL-BEING

Inpatient and outpatient therapy



www.gardenbalance.it



Two separate back and seat/leg cushions with removable seat/leg cushion.



Multi-adjustable ergonomic armrests.

Wheels with braking system activated by bilateral levers.

Base protection covers made of ABS, with a smooth surface and easy to sanitise.




Multi-adjustable ergonomic armrests

Ergonomic armrests equipped with a special manual mechanism with multiple adjustment for positioning in different directions.

In addition, the support placed on both sides of the armrests (accessory code AC25P) allows the hand control to be conveniently stored when not in use.



Removable

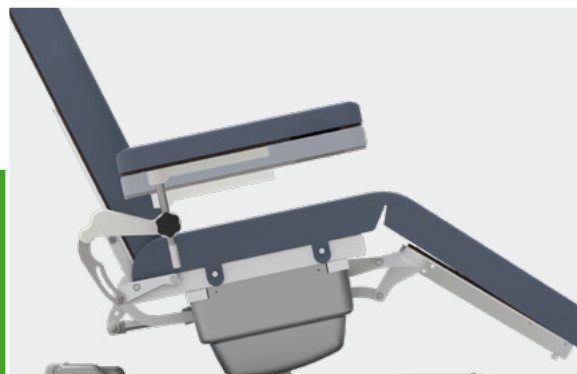

Adjustable height



Solidal to the Backrest



Rotating



Tiltable



Overturned at 90°




Movements

Synchronised back and leg sections

The movements of the synchronised back and leg sections are controlled by dedicated buttons on the handset.



Activated by
Handset

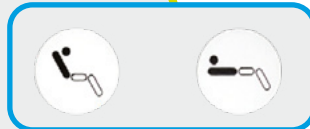

Movements simultaneous
of the back and leg sections




Movements

Trendelenburg


The Trendelenburg position can be obtained using the dedicated buttons on the handset.

Activated by
Handset

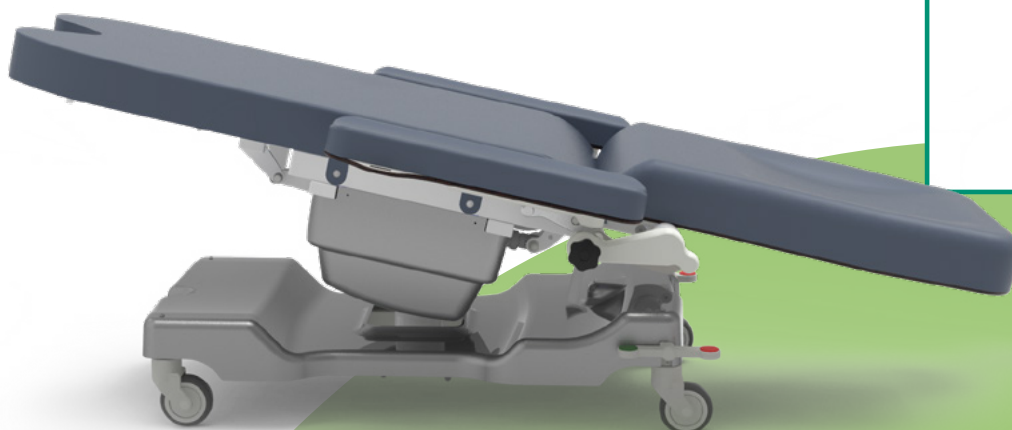

Shock position

Another example of movement

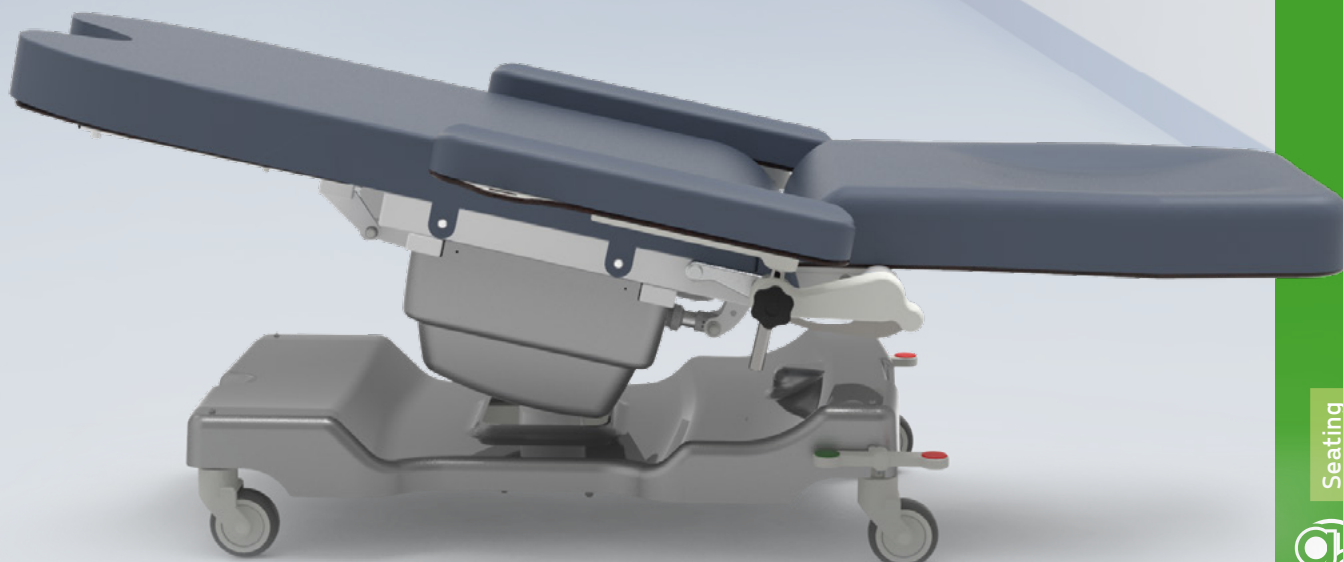



Other features of the medical device

- Painted steel frame in thermosetting powder
- Lying surface divided in 3 sections (backrest, seat and leg).
- Simultaneous movement of back and leg section through 1 motor.
- Trendelenburg position (shock position as alternative).
- Armrest's movements synchronized with the backrest. The armrests are completely removable and can be rotate, tilted and overturned at 90°; (please refer to the accessories catalogue for other version of armrests).
- No. 4 wheels ø 100 mm antistatic, sturdy, stainless and swivable.
- Simultaneous brake on rear wheels in chromed steel - Equipment code AC7P1G.
- Chair's movements are controlled by handset with a 3,5 mtl coiled cable. The handset is user friendly and has visual pictograms useful for partially sighted people. Handset is equipped by 10 buttons with 2 memories.
For patient's safety purpose the handset is provided with an inhibition key - Equipment code AC34P. Insulation IP66.
- Bilateral handset support - Equipment code AC25P.
- ABS base cover - Equipment code AC83.
- ABS seat cover - Equipment code AC84.
- Easily removable soft foam independent cushions. Upholstery in flame retardant leatherette UNI 9175/87 e 9175FA-1/94 class 11M - Equipment code AC87/1.

Do you want to change the colour of the upholstery and frame?

Frame or click the QR codes and discover how to personalise your device by consulting our colour cards*.

*If no preference is expressed, the device will be delivered in the available colour.

Upholstery
Leatherette

Coating
Silicon

Frame



Do you want to improve the device?

Scan the QR Code to find out all the available accessories listed on our catalogue.



Technical sheet

CND/EMDN code	V08020101
GMDN Code	38447
N. progressive/R	2433992/R
UDI-DI di Base	805771740POLPE
ID Product	GB0165.SP Painted steel
Inteded use	Multifunctional chair
Manufacturer	GARDHEN BILANCE SRL
Overall weight	Standard 70 kg / Bariatric 100 kg
Wheels	No.4 Ø 100 mm swivable and stainless with braking system
Power	100 - 240 Vac 50/60 Hz
Electric motors	Bassa tensione (24 V)
Max electrical input	4A
Insulation	IP44
Max load	Standard 200 kg / Bariatric 300 kg
Safety working load	Standard 235 kg / Bariatric 335 kg

Main cord



Type F
Schuko



Type L
16A



Type I
Australia





Type B
Canada and USA



Type G
UK plug



Dimensions

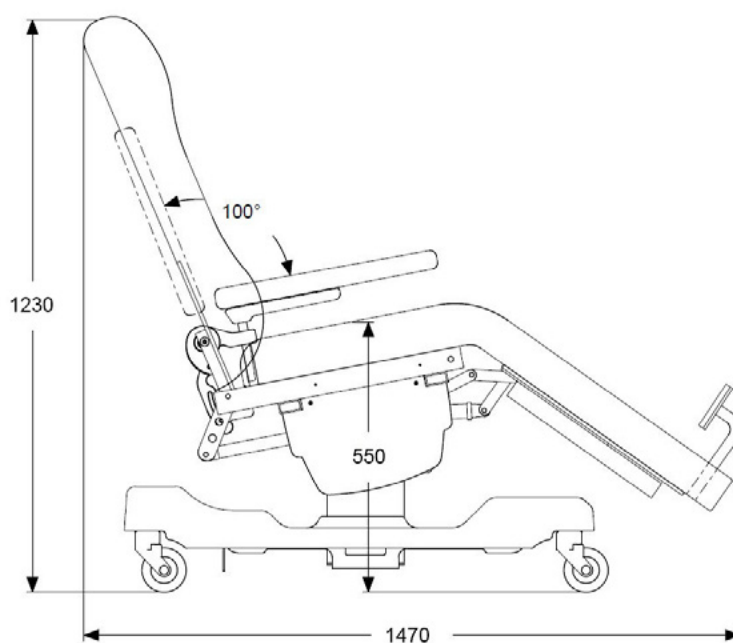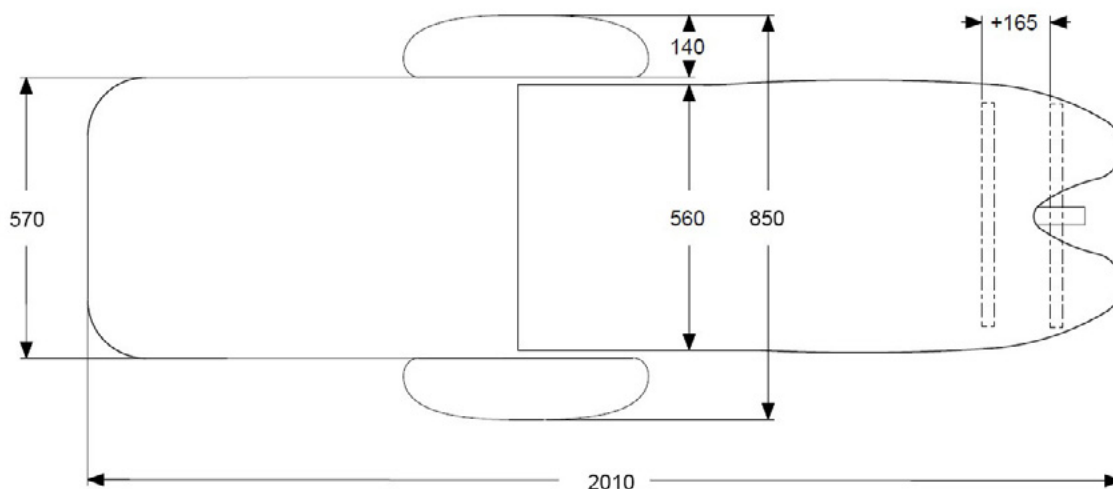

www.gardhenbalance.it



Gardhen Balance S.r.l.

Via G. Luraghi c/o Consorzio Il Sole - lotto S - 80038
Pomigliano d'Arco (NA) - Italy

Tel. +39.081.8692160 pbx - Fax +39.081.8692460

E-mail: info@gardhenbalance.it

PEC: gardhenbalance@pec.it

Pictures contained in this brochure is for illustrative purpose only. The manufacturer reserves the right to modify the projects without notice. For accessories identification please refer to the accessory catalogue or commercial offer. When data are different please refer to commercial offer.

Dep. EN_88 rev. 16/10/2024