

Cloud - Electric Stretcher

VARIABLE HEIGHT



STRETCHERS



Frequent Usage

Patients Transport
Emergency room

it is a product



FUNCTIONAL & SAFE


Maximum comfort
at the service of safety



Partner and Distributor
Pam Mobility



Cloud - Electric stretcher

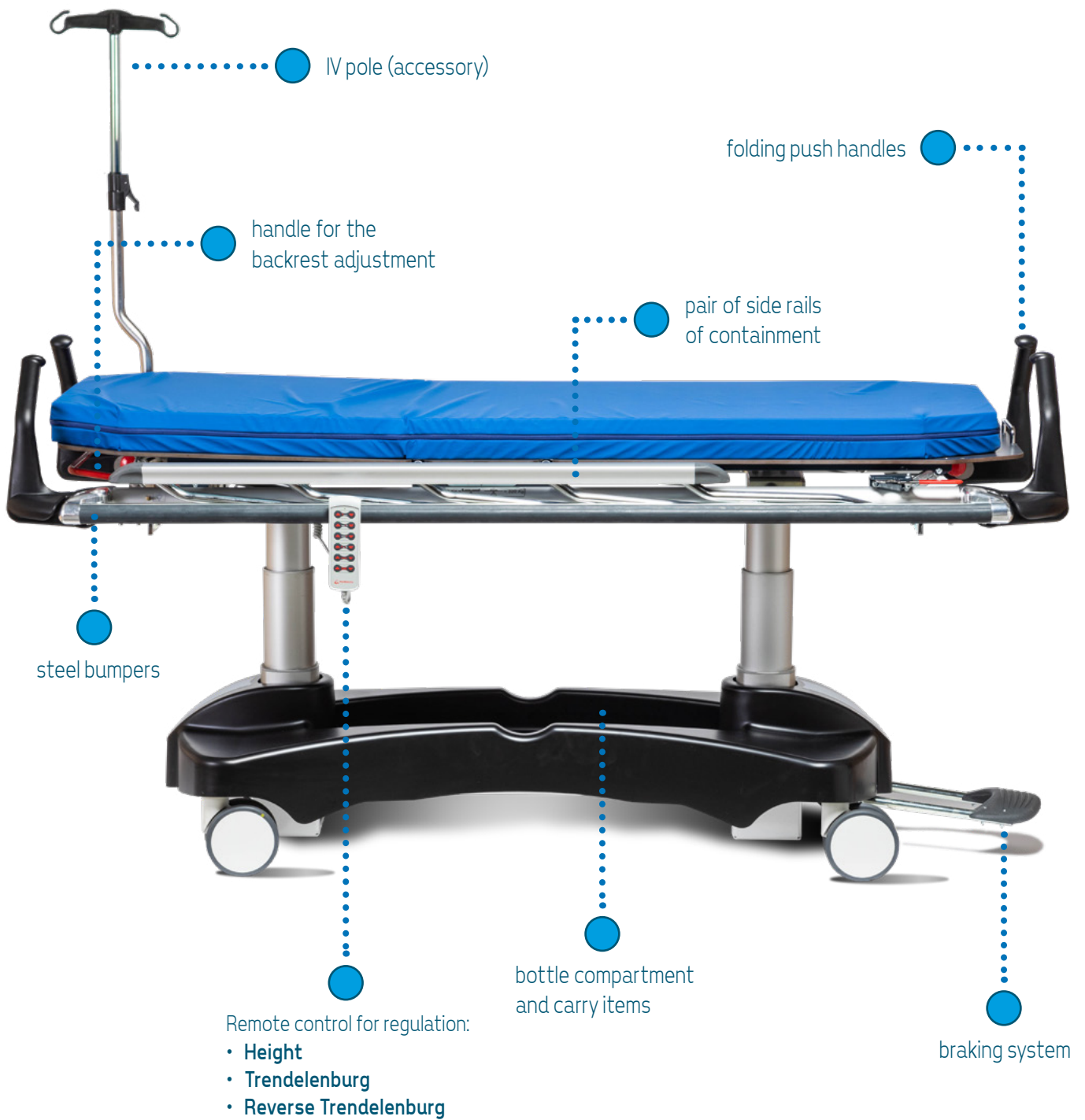



Configurations of the sections of the lying plane



Model 7285RT70 - 2 sections with increased top

- Backrest section adjustable by gas pistons.
- Fixed basin section.
- Adjustable leg section with 5-position telescopic lifting joint.



Model 7285RT71 - 4 sections with increased top

- Backrest section adjustable by gas pistons
- Fixed basin section.
- Femoral section adjustable by gas pistons using the handle placed at the end of the foot section.
- Adjustable leg section with 5-position telescopic lifting joint.

Standard equipment common to all models

- Steel mesh frame with integrated preparation of multifunctional supports at the corners, which allow the insertion of accessories.
- All steel parts are painted with epoxy powders, easy to sanitize and resistant to corrosion and chemical agents.
- HPL radiolucent lying surface with mattress stops.
- Folding steel handles placed on the foot end.
- Push handle in steel placed on the head side.
- Steel bumpers placed at the ends of the frame and on the perimeter.
- Variable height electrically adjustable via push-button panel, bringing the device from a minimum height of 500 mm up to a maximum height of 900 mm.
- Motorized trendelenburg and reverse trendelenburg ($- 18^\circ / + 18^\circ$) adjustable via push-button panel.
- Nr. 4 wheels \varnothing 200 mm antistatic, non-marking, revolving on ball bearings.
- 5th \varnothing 150 mm directional wheel activated with pedal.
- Braking system with 3-position foot pedal:
 - Free wheels;
 - Braked wheels;
 - 5th directional wheel.
- Pair of lateral containment sides made of compass-shaped extruded aluminum alloy (collapsible), easy to sanitize.
- Accessories holder bar.
- Protective shaped technopolymer cover for all mechanisms and flat surfaces.
- Compartment located on the cylinder holder cover and/or container compartment for inserting other objects as needed.
- Static mattress, with removable cover made of open cell expanded polyurethane foam with





a height of 6 cm

- Mattress cover made with a special bi-elastic fabric (washable, hypoallergenic and fireproof class 1M with antibacterial and antifungal properties and resistant to sanitization with common disinfectants used in hospitals and outpatient clinics).
- The sanitization method of the device offered is compatible with (non-dedicated) chlorinated products capable of developing an amount of available chlorine (approximate to active chlorine) equal to 500-1000 ppm or hydroalcoholic solutions based on chlorhexidine gluconate at the 0.5% or hydrogen peroxide-based products (e.g. Incidin OXY FOAM flc. 750 ml for surfaces).



**FUNCTIONAL &
SAFE**



Movements:

The **Trendelenburg** and **reverse Trendelenburg** movement is motorized and adjustable with push button on panel.

TRENDELENBURG ■



REVERSE
TRENDELENBURG ■





Movement:

The **variable height of the stretcher** is adjustable by means of a pedal placed on both sides, bringing the device from a minimum of 500 mm up to a maximum of 900 mm.





Configuration

Model	large lying sections	
	7285RT70	7285RT71
width surface (cm)	70	70
length surface (cm)	195	195
overall width (cm)	87	87
overall length (cm)	207,5	207,5
Number of sections	3	4
backrest inclination (degrees)	80°	80°
femoral inclination (degrees)	-	21°
leg inclination (degrees)	-	20°
trendelenburg/reverse trendelenburg (degrees)	±18°	±18°
max load capacity (kg)	260	260
safety working load (kg)	320	320
minimum height (mm)	500	500
maximum height (mm)	900	900




Do you want to implement the device?

Frame the QR code and discover all the available accessories by consulting the catalogue.



Dimensions


Model 7285RT70 - 2 sections with large lying sections

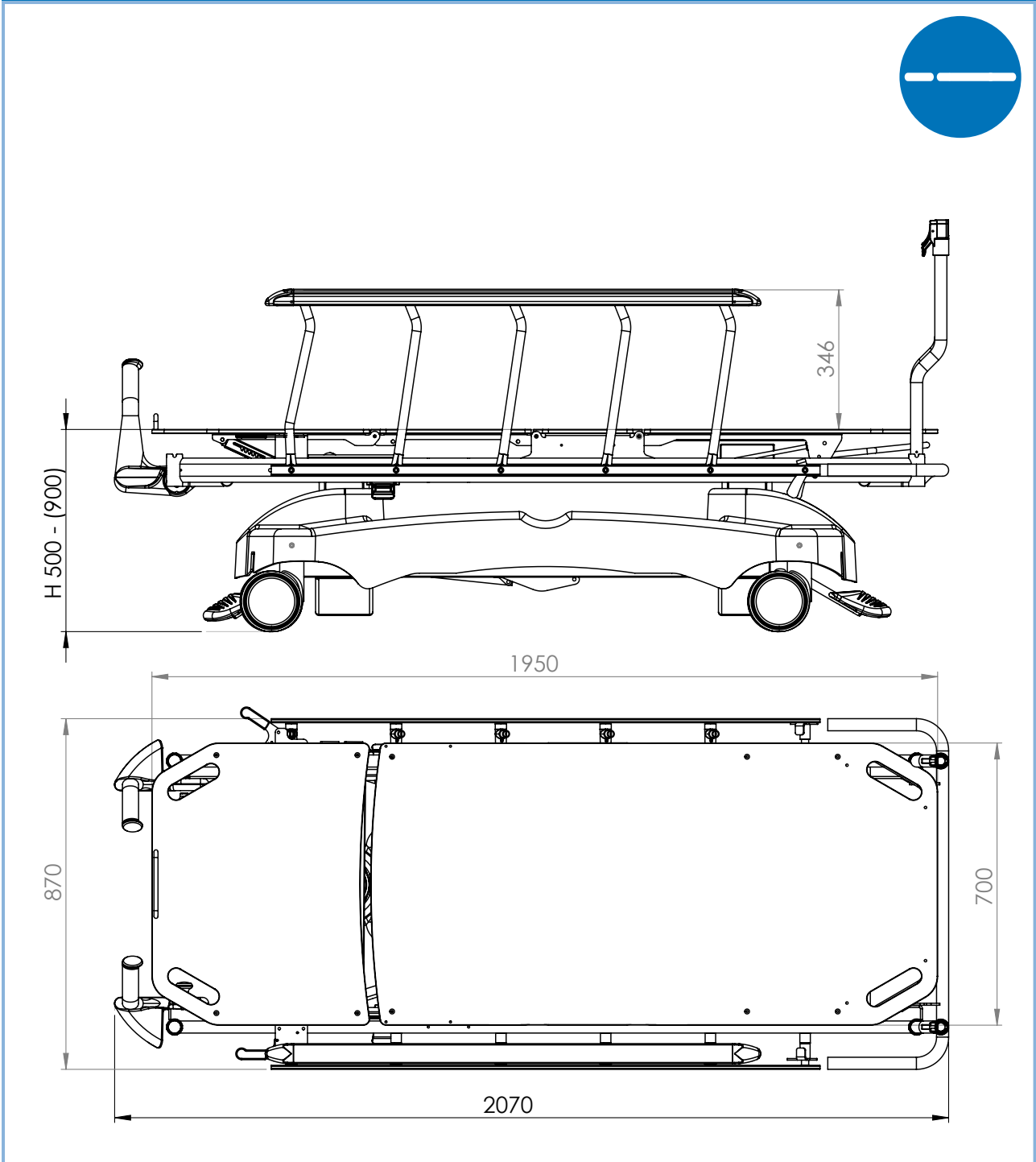



Dimensions

Model 7285RT71 - 4 sections with large lying sections

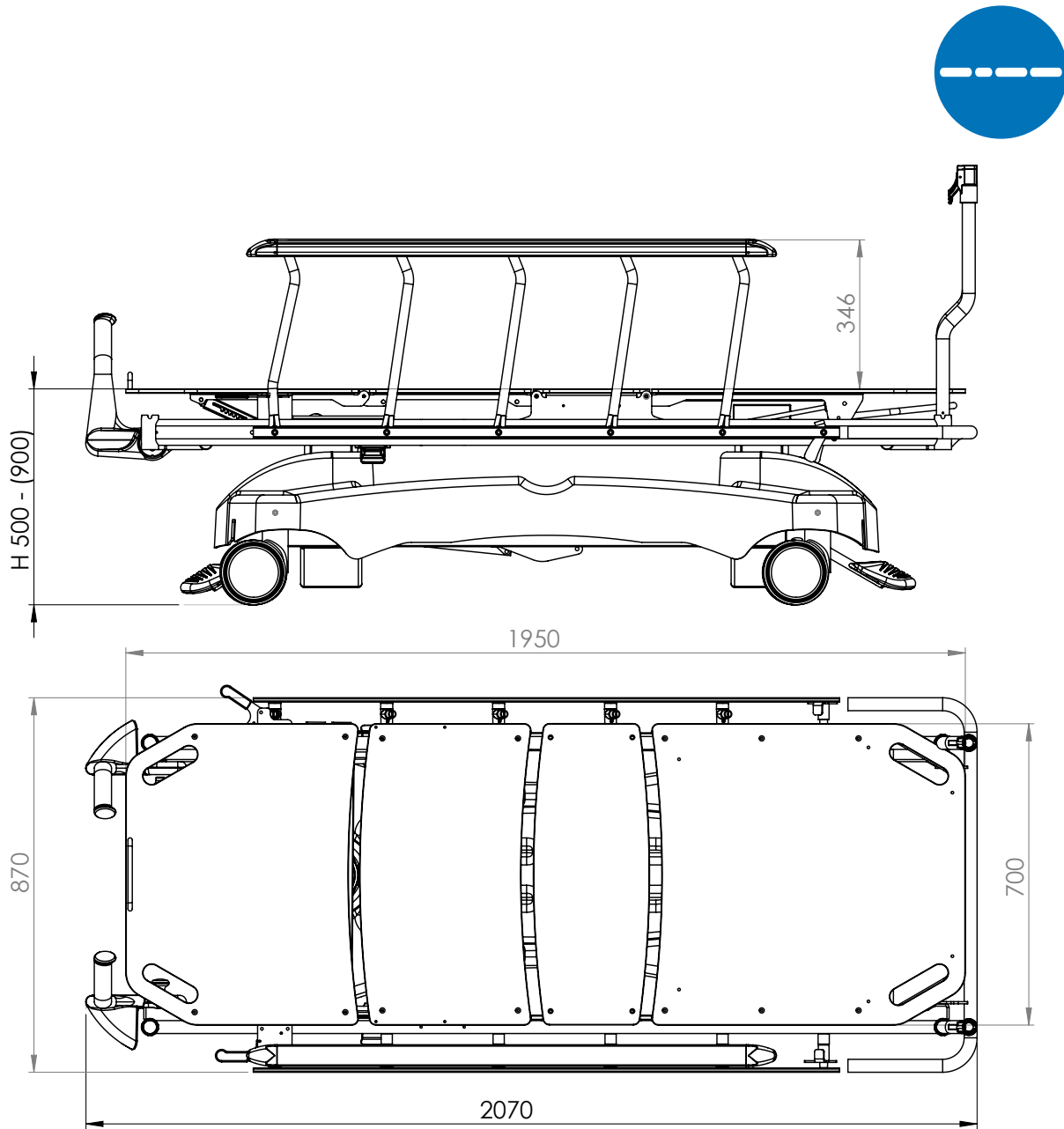



Technical sheet

ID code	7285RT70 - 2 sections with large lying sections
ID code	7285RT71 - 4 sections with large lying sections
Intended use	Stretcher for patient's transport

It is a  **PamMobility** product

Gardhen Bilance is a partner and exclusive distributor for:

THE INTERNATIONAL MARKET EXCEPT IN: Germany, Sweden, Finland, Austria, Belgium, Switzerland, the Netherlands, UK.