

# Cloud - Electric Stretcher

VARIABLE HEIGHT



STRETCHERS



Frequent Usage

Patients Transport  
Emergency room

it is a product

 **PamMobility**

# FUNCTIONAL & SAFE


Maximum comfort  
at the service of safety



Partner and Distributor  
Pam Mobility



# Cloud - Electric stretcher

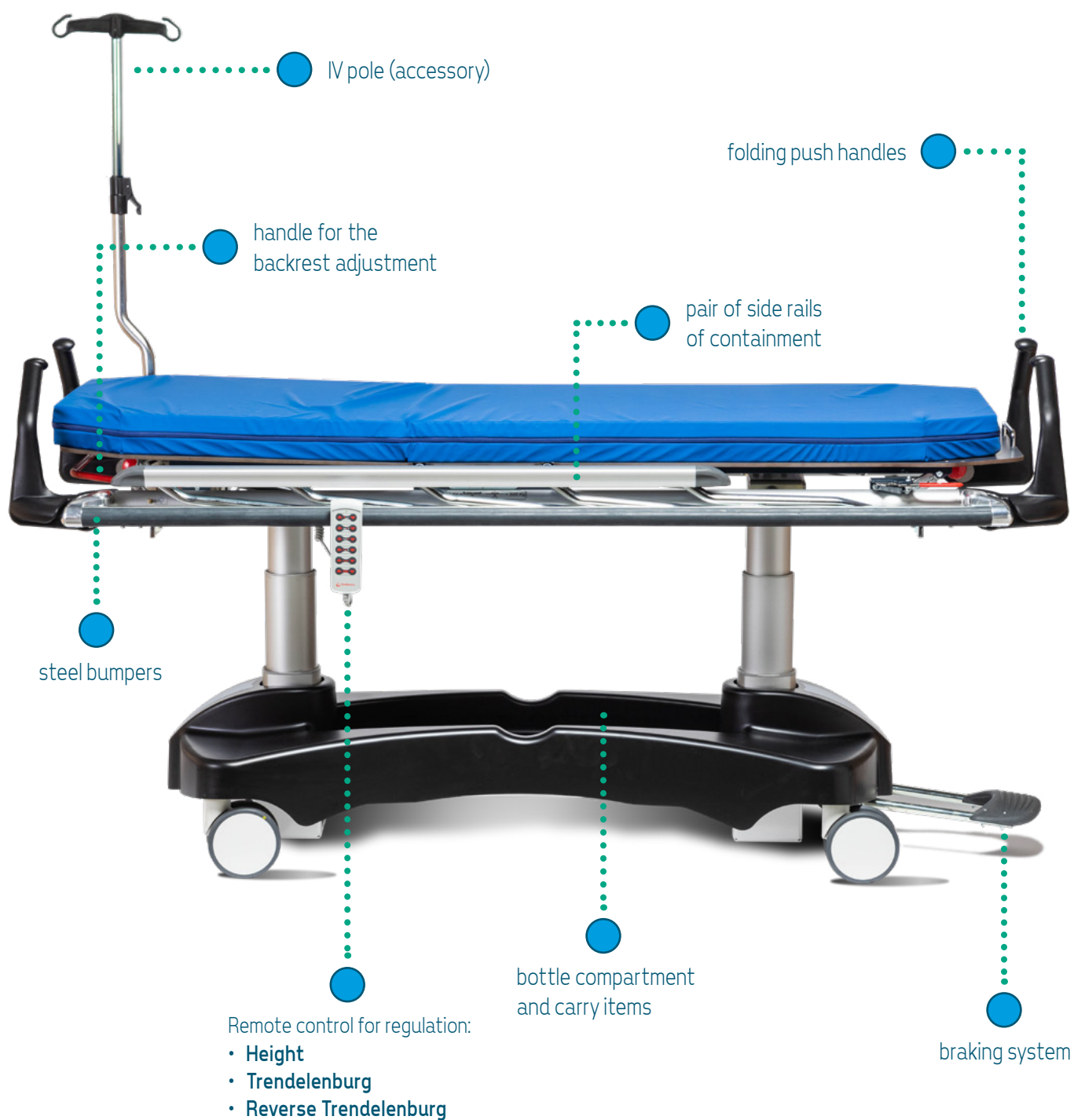



## Configurations of the sections of the lying plane



### **Model 7285RT70 - 2 sections with increased top**

- Backrest section adjustable by gas pistons.
- Fixed basin section.
- Adjustable leg section with 5-position telescopic lifting joint.



### **Model 7285RT71 - 4 sections with increased top**

- Backrest section adjustable by gas pistons
- Fixed basin section.
- Femoral section adjustable by gas pistons using the handle placed at the end of the foot section.
- Adjustable leg section with 5-position telescopic lifting joint.

## Standard equipment common to all models

- Steel mesh frame with integrated preparation of multifunctional supports at the corners, which allow the insertion of accessories.
- All steel parts are painted with epoxy powders, easy to sanitize and resistant to corrosion and chemical agents.
- HPL radiolucent lying surface with mattress stops.
- Folding steel handles placed on the foot end.
- Push handle in steel placed on the head side.
- Steel bumpers placed at the ends of the frame and on the perimeter.
- Variable height electrically adjustable via push-button panel, bringing the device from a minimum height of 500 mm up to a maximum height of 900 mm.
- Motorized trendelenburg and reverse trendelenburg adjustable via push-button panel.
- Nr. 4 wheels ø 200 mm antistatic, non-marking, revolving on ball bearings.
- 5<sup>th</sup> ø 150 mm directional wheel activated with pedal.
- Braking system with 3-position foot pedal:
  - Free wheels;
  - Braked wheels;
  - 5<sup>th</sup> directional wheel.
- Pair of lateral containment sides made of compass-shaped extruded aluminum alloy (collapsible), easy to sanitize.
- Accessories holder bar.
- Protective shaped technopolymer cover for all mechanisms and flat surfaces.
- Compartment located on the cylinder holder cover and/or container compartment for inserting other objects as needed.
- Static mattress, with removable cover made of open cell expanded polyurethane foam with a height of 6 cm



- Mattress cover made with a special bi-elastic fabric (washable, hypoallergenic and fireproof class 11M with antibacterial and antifungal properties and resistant to sanitization with common disinfectants used in hospitals and outpatient clinics).
- The sanitization method of the device offered is compatible with (non-dedicated) chlorinated products capable of developing an amount of available chlorine (approximate to active chlorine) equal to 500-1000 ppm or hydroalcoholic solutions based on chlorhexidine gluconate at the 0.5% or hydrogen peroxide-based products (e.g. Incidin OXY FOAM flc. 750 ml for surfaces).



**FUNCTIONAL &  
SAFE**



## Movements:

The **Trendelenburg** and **reverse Trendelenburg** movement is motorized and adjustable with push button on panel.

TRENDELENBURG ■



REVERSE  
TRENDELENBURG ■







### Movement:

The **variable height of the stretcher** is adjustable by means of a pedal placed on both sides, bringing the device from a minimum of 500 mm up to a maximum of 900 mm.





## Configuration

Model	large lying sections	
	7285RT70	7285RT71
width surface (cm)	70	70
length surface (cm)	195	195
overall width (cm)	87	87
overall length (cm)	207,5	207,5
Number of sections	3	4
max load capacity (kg)	260	260
safety working load (kg)	320	320
minimum height (mm)	500	500
maximum height (mm)	900	900




### Do you want to implement the device?

Frame the QR code and discover all the available accessories by consulting the catalogue.



## Dimensions

### Model 7285RT70 - 2 sections with large lying sections

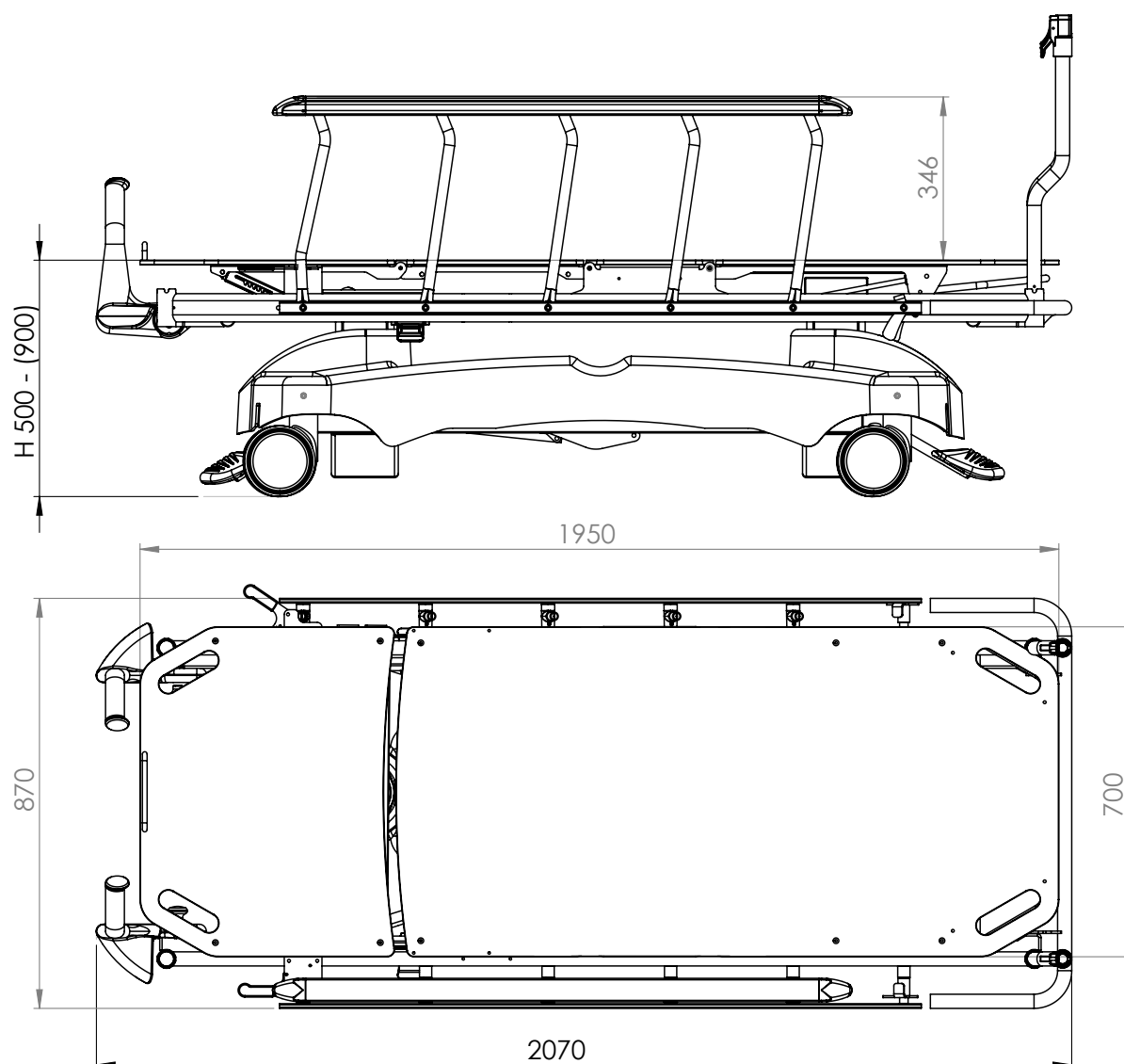






## Dimensions

### Model 7285RT71 - 4 sections with large lying sections

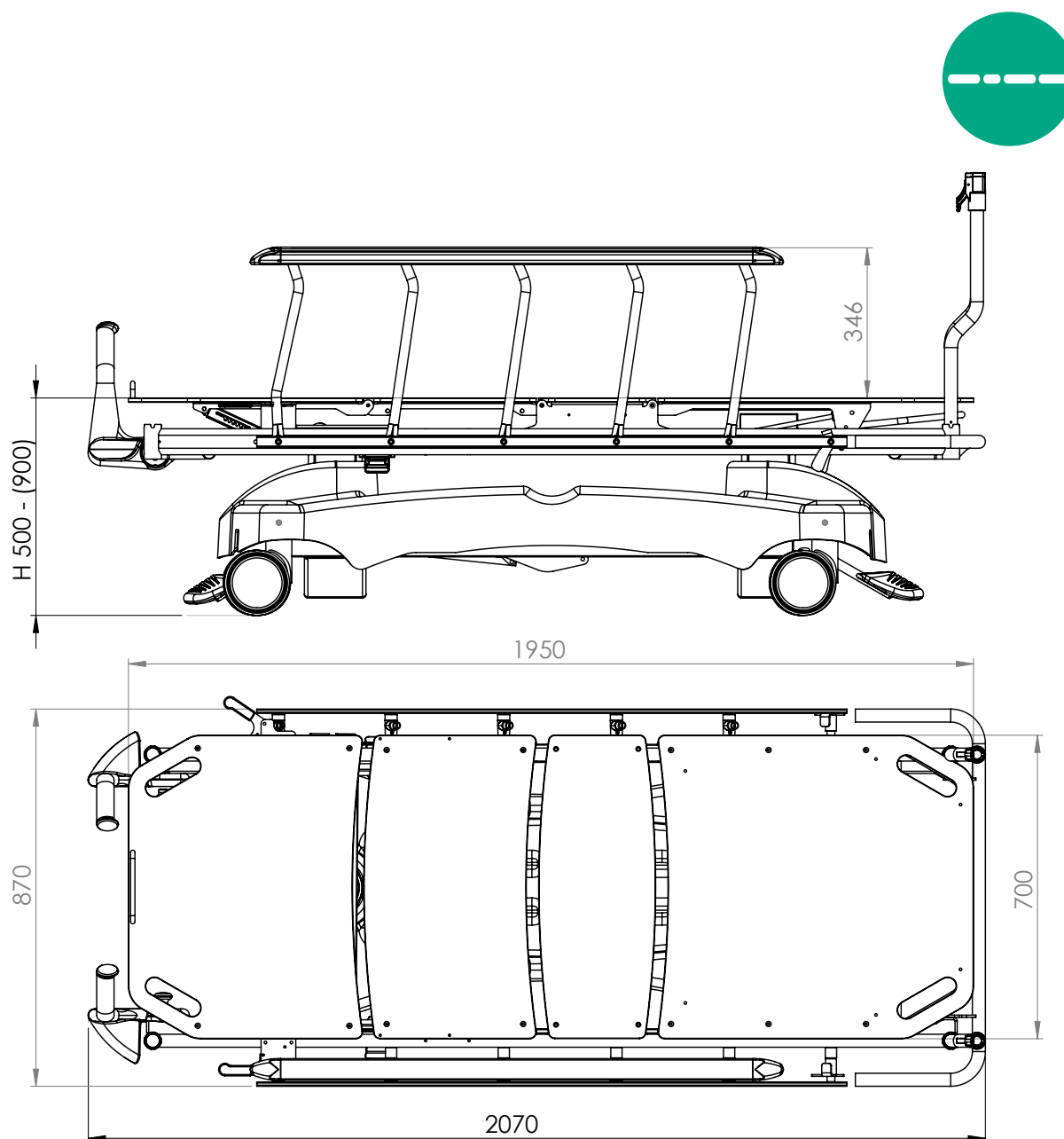



### Technical sheet

ID code	7285RT70 - 2 sections with large lying sections
ID code	7285RT71 - 4 sections with large lying sections
Intended use	Stretcher for patient's transport

It is a  **PamMobility** product

**Gardhen Balance is a partner and exclusive distributor for:**

THE INTERNATIONAL MARKET EXCEPT IN: Germany, Sweden, Finland, Austria, Belgium, Switzerland, the Netherlands, UK.