

Flex

Gynecological/urological examination table



COUCHES

Frequent use

Gynecology
Urology
Outpatient clinic
Diagnostic Ultrasound
Endoscopic Examinations
Medical Centre
Emergency Room



**VERSATILE &
SAFE**

multifunctional device
for specialized activities



www.gardhenbalance.it



VERSATILE & SAFE

The **examination couch Flex** is extremely versatile and can be used for various gynecological/urological examinations. With its three sections and ergonomic thigh supports, it allows for optimal patient posture while ensuring greater comfort during the examination or other medical procedures.

ROLL HOLDER

Made of steel placed on the back of the backrest.

LEG BRACE

Shaped manually adjustable.

WHEELS

Which allow the device to be easily moved from place to place.

FEET

Which allow the device to be placed stably on all pavements.

BASIN

Made of stainless steel for liquid collection.





Standard equipment

- Electric variable height gynecological/urological examination table.
- Reinforced frame made of steel painted with thermosetting epoxy powder enriched with antibacterial, color customization upon request.
- Table size 1850 mm W - 650 mm D.
- Top divided into 3 sections and 2 joints. Leg section with gas-spring adjustable recess. Leg section foldable even with patient positioned on the bed.
- Manually adjustable back section with piston (gas spring), movement operable with the help of lever placed on both sides, to have a functional tilt to the type of visit/therapy.
- Electrically adjustable variable height, minimum height 580 mm, maximum height 830 mm, with total range of 250 mm.
- Laying surface (80 mm in height) composed of 3 cushions made of differentiated density polyurethane. All the materials making up the padding are honeycomb structure to ensure high breathability.
- Upholstery made of class 11M leatherette, a material that can be sanitized with commonly used detergents and medical surgical aids, with antibacterial, antifungal, stain-resistant properties. Color customization upon request - Equipment code AC87/1.
- Quick release of the backrest, to facilitate CPR maneuver - Equipment code AC8L.
- Contoured, 3D adjustable, soft and comfortable leg brace for patients- Equipment code AC110/1.
- Removable stainless steel basin for collecting liquids - Equipment code AC155.
- 2 wheels and 2 non-slip height-adjustable feet.
- Perimeter 2-button push-button panel for adjusting the height of the top.
- Painted steel roll holder placed on the back of the backrest - Equipment code AC11.





ELECTRICALLY ADJUSTABLE VARIABLE HEIGHT




Do you want to implement the device?

Frame or click the QR code and discover all available accessories by checking out the **accessories tab**.



Dimensional

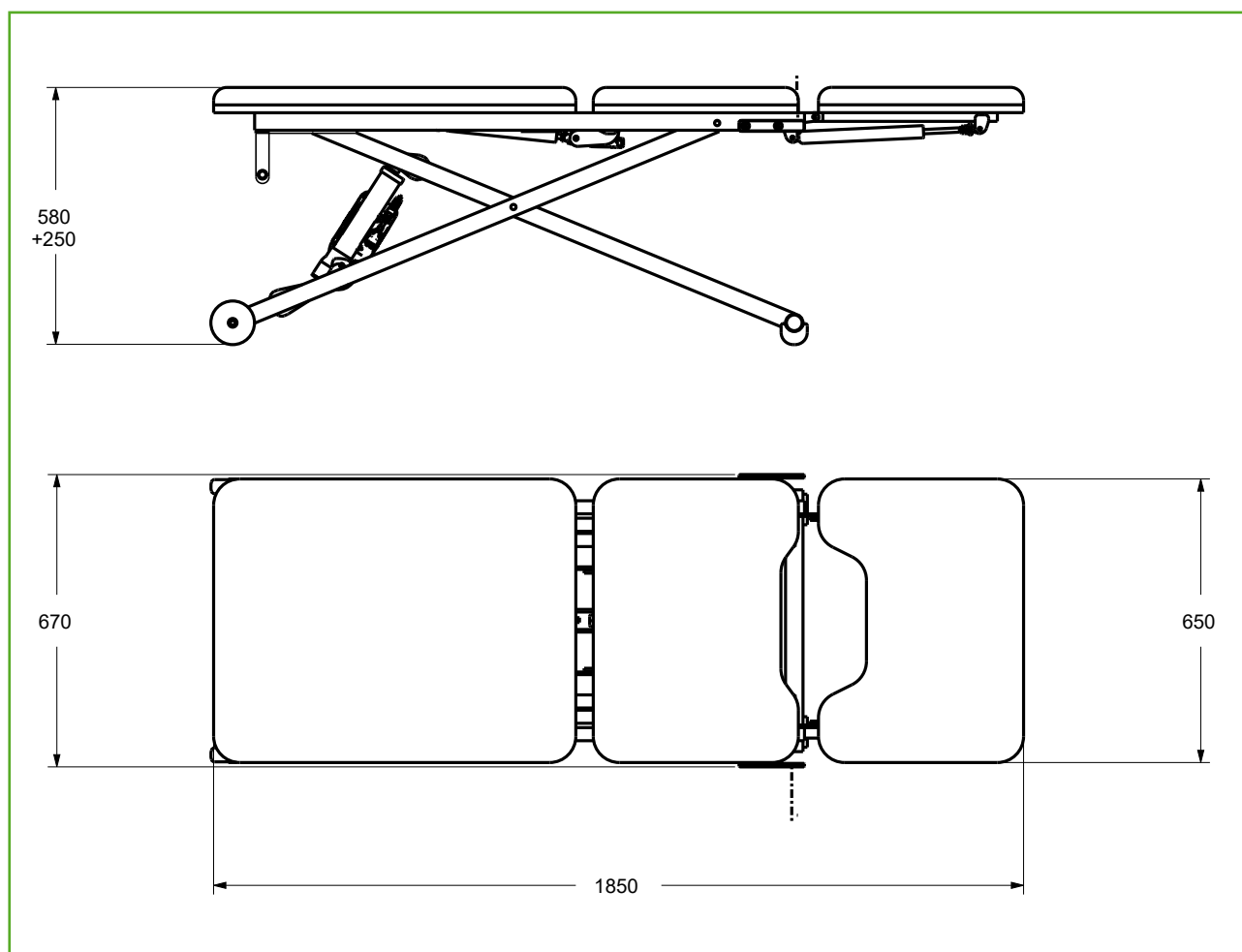

Do you want to change the color of the frame?

Frame or click the QR code and find out how to personalize the device by consulting our color card.



Do you want to change the color of the coating?

Frame or click the QR code and find out how to personalize the device by consulting our color card.



Data sheet

CND - EMDN Code	V080603
GMDN Code	34870
No. progressive/R	2402702/R
UDI-DI di base	805771740LVISKK
ID Product	GB0400.SP
Intended Use	Gynecological/urological crib
Overall weight	70 kg
Type of control	Push button panel
Power supply	Multi-voltage 100-240 Vac 50/60 Hz
Power supply with batteries	24v - 4.5A
Electric motors	Low voltage (24 V)
Max electrical input	4A
Insulation	IP44
Max load capacity	200 kg
Safety working load	235 kg

Main cord



Type F - Schuko



Type L - 16A

Gardhen Balance S.r.l.

Via G. Luraghi c/o Consorzio Il Sole - lotto S - 80038 Pomigliano d'Arco (NA) - Italy

Tel. +39.081.8692160 pbx - Fax +39.081.8692460

E-mail: info@gardhenbalance.it - PEC: gardhenbalance@pec.it

Visit our website www.gardhenbalance.it / Follow us on     

Images reproduced in this sheet are for reference only and dimensional values are nominal. It is the manufacturer's right to make changes to the designs without prior notice. Refer to the accessories catalog and/or commercial proposals for the identification of standard equipment and the selection of any additional components. In case of discrepant data, the contents of the commercial proposals shall prevail. For regulatory references, please refer to the user manual and declaration of conformity.
Dep. 336 rev. 30/04/2024