

Physio

Physiotherapy and massage table
electric variable height



COUCHES

Frequent usage

Physiotherapy
Rehabilitation
Diagnostics
Medical Clinic

**DOUBLE
BACKREST**

**VERSION WITH
OPTIONAL
COMPONENTS**

- Roll holder
- Facial Hole

**VERSATILE &
SAFE**

professional support for medical treatment



www.gardhenbalance.it



VERSATILE & SAFE

The **physiotherapy and massage table** is extremely versatile thanks to the adjustment of the two backrests.

- **Low backrest**, adjustable positive and negative, can be equipped with face hole upon request.
- **High backrest** for people to assume sitting positions.

Various configuration options ensure maximum patient comfort during examination or other medical procedures.

HIGH BACKREST

Adjustable to make people assume the sitting position

LOW BACKREST

Adjustable to optimize patient posture.

DOUBLE BACKREST



Standard equipment

- Top divided into 3 sections - 2 joints.
- Reinforced frame made of steel painted with thermosetting epoxy powder enriched with antibacterial, color customization upon request.
- Low back/leg adjustable with piston (gas spring), the movement is manually operated with the help of lever placed on both sides.
- High back/leg adjustable with piston (gas spring), the movement is manually operated with the help of lever placed on both sides.
- Trendelenburg obtainable with the help of levers placed on either side.
- Height adjustment obtained by pantograph compasses operated electrically via push-button panel. Minimum height 510 mm, maximum height 910 mm, with 400 mm range.
- Perimeter handset equipped with 3.5 m coiled cable, The handset is easy to use as it is made with large raised pictograms and is also suitable for visually impaired patients; IP66 insulation.
- 4 swivel casters with brakes (Ø 125 mm), non-marking, low rolling resistance, high quietness, sliding on bearings and high capacity. Each wheel is equipped with special brake that can be operated by foot lever.
- Padding made of high-density foam resin.
- Upholstery made of class 11M leatherette, a material that can be sanitized with commonly used detergents and medical surgical aids, with antibacterial, antifungal, stain-resistant properties. Color customization upon request - Equipment code AC87/1.

Do you want to change the color of the upholstery and frame?

Frame or click the QR codes and find out how to personalize the device by consulting our color cards*.

*If no preference is expressed, the device will be supplied with the available color.

**Leatherette
Upholstery**

**Siliconic
Upholstery**

Frame




Do you want to implement the device?

Frame the QR code and discover all the optional components available by browsing the catalog.



Movements

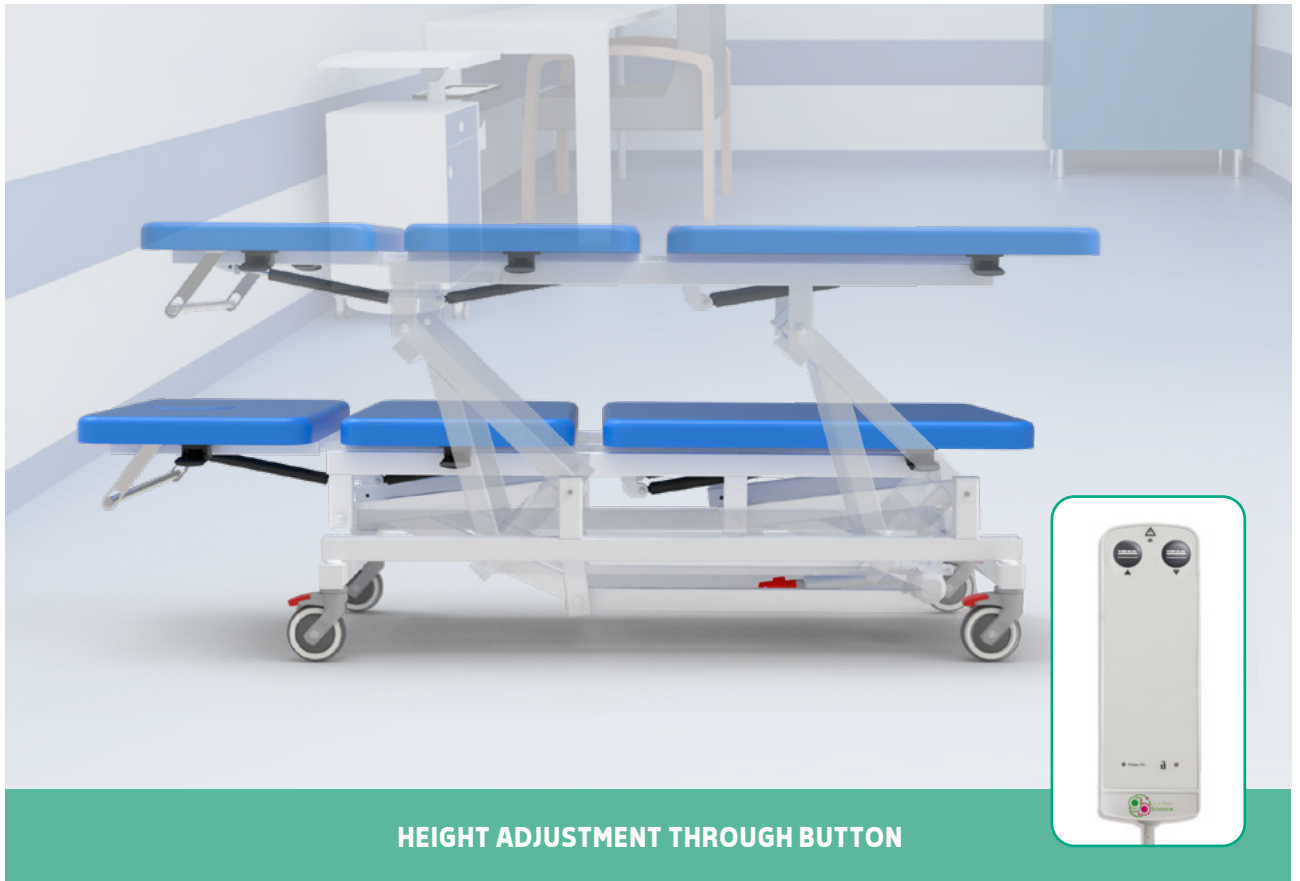

HEIGHT ADJUSTMENT THROUGH BUTTON



LOW BACK ADJUSTMENT WITH PISTON (GAS SPRING)




HIGH BACK ADJUSTMENT WITH PISTON (GAS SPRING)



TRENDELENBURG OBTAINABLE WITH THE AID OF LEVERS PLACED ON EITHER SIDE



Dimensions

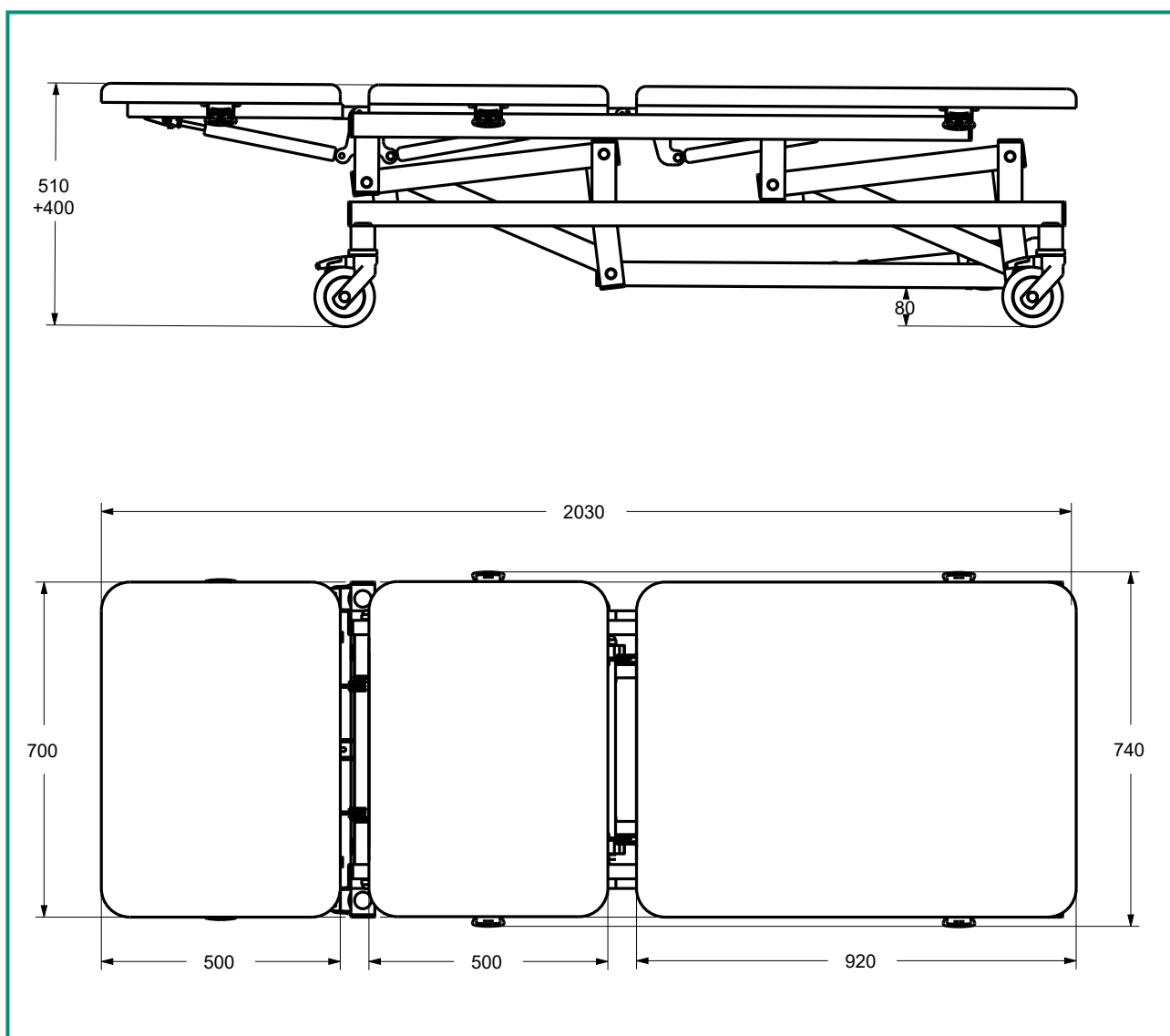

Main cord



Type F
Schuko



Type L
16A



Type I
Australia



Type B
Canada
and USA



Type G
UK plug



Data sheet

CND/EMDN	V080603
GMDN	34870
UDI-DI di Base	805771740LVISKK
No. progressive/R	2293770/R
Product Code	GB0389.SP
Intended Use	Physiotherapy/massage table
Manufacturer	GARDHEN BILANCE SRL
Total weight	70 kg
Type of control	Push button panel
Ruote	No. 4 Ø 125 mm, stainless swivel with braking system
Power supply	Multi-voltage 100-240 Vac 50/60 Hz
Power supply with batteries	24v - 4.5A
Electric motors	Low voltage (24 V)
Max absorption	4A
Insulation	IP44
Max load	180 kg
Safety working load	215 kg

Gardhen bilance

Via G. Luraghi c/o Consorzio Il Sole - lotto S - 80038 Pomigliano d'Arco (NA) - Italy

Tel. +39.081.8692160 pbx - Fax +39.081.8692460

E-mail: info@gardhenbilance.it - PEC: gardhenbilance@pec.it

Visit our website **www.gardhenbilance.it** / Follow us on     

Images reproduced in this sheet are for reference only and dimensional values are nominal. It is the manufacturer's right to make changes to the designs without prior notice. Refer to the optional components catalog and/or commercial proposals for identification of standard equipment and selection of any additional components. In case of discrepant data, the contents of the commercial proposals shall prevail. For regulatory references, please refer to the user manual and declaration of conformity.

Dep. EN_316 rev 29/09/2024