

Stephen H Eco 2.0

Multifunctional armchair



SEATING

Frequent Usage

Chemotherapy
Infusion therapy
Day Hospital
Blood Donation
Blood samples
Emergency room
Dialysis
Patient transport
Dermatology/
Cosmetics
Multifunctional



High

COMFORT

Inpatient and outpatient therapy




www.gardhenbalance.it



Standard equipment

Ergonomic armrests at multiple adjustment equipped with support for the hand control

(Equipment code AC25P) placed bilaterally on the armrests.

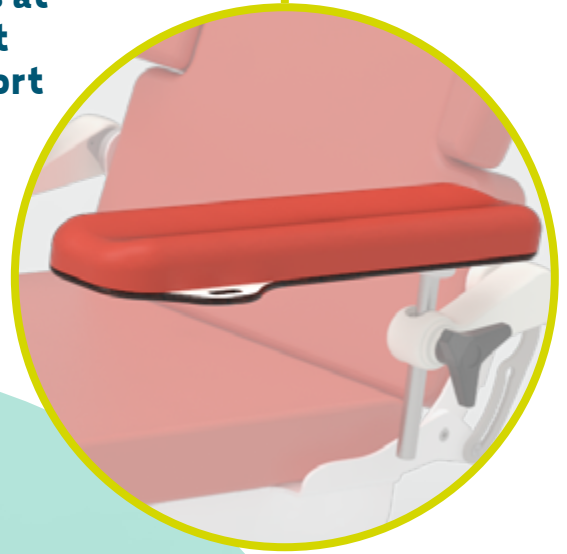

Handle made of stainless steel to facilitate carrying activities. Equipment code AC56.





Standard equipment



Supportive of the Backrest



Swivels



Height adjustable



Removable



Tiltable



Tipplers



Movements

Handset that controls the movements provided.

Movements simultaneous of the backrest and leg sections

Trendelenburg





Characteristics of the medical device

- Steel frame painted with thermosetting epoxy powder enriched with thermosetting powder antibacterial.
- Anti-tip frame ensuring stability under all operating conditions.
- Top divided into 3 sections: backrest sections, seat section and legs section.
- Simultaneous movement of back and leg sections by 1 motor.
- Trendelenburg position (alternatively shock position).
- Armrests integral to the backrest movement, with swivel function, tiltable up/down, 90° tilt and height adjustment; (see optional components catalog for other armrest type).
- No. 4 durable, antistatic, stainless, swivel casters Ø 100 mm.
- Simultaneous brake galvanized steel rear wheels - Equipment code AC7P1G.
- All chair movements are controlled by a perimeter handset equipped with a 3.5 m cable. The handset is easy to use as it is made with large raised pictograms and is also suitable for visually impaired patients. The handset has 10 keys with 2 memories. For patient safety, the handset is equipped with a key for control inhibition - Equipment code AC34P. IP66 insulation.
- bilateral support on the armrests for the handset - Equipment code AC25P.
- Stainless steel carrying handle-Equipment code AC56.
- Two separate back and seat/leg cushions, made of functional soft foam for: excellent posture, ergonomics and preventive effect of decubitus injuries.
- Upholstery made of class 11M leatherette, a material that can be sanitized with commonly used detergents and medical surgical aids, with antibacterial, antifungal, stain-resistant properties. Color customization upon request - Equipment code AC87/1.
- Fully removable liner to enable effective and functional sanitization and to reduce cross-infection.
- Device that can be sanitized with commonly used medical surgical supplies and detergents.

Do you want to change the color of the upholstery and frame?

Frame or click the QR codes and find out how to customize your device by browsing our color cards*.

*If no preference is expressed, the device will be supplied with the available color.

Coating
Similpelle

Coating
Silicon

Frame





Do you want to implement the device?

Frame the QR code and discover all the optional components available by browsing the catalog.



Size table

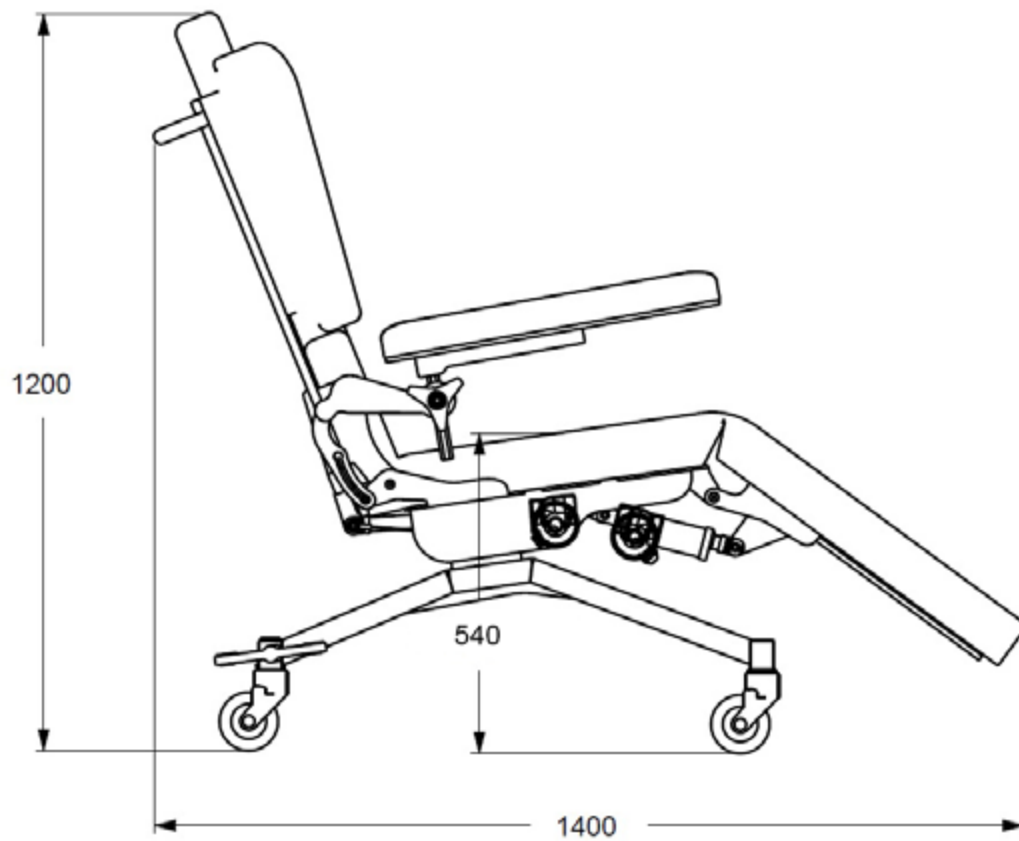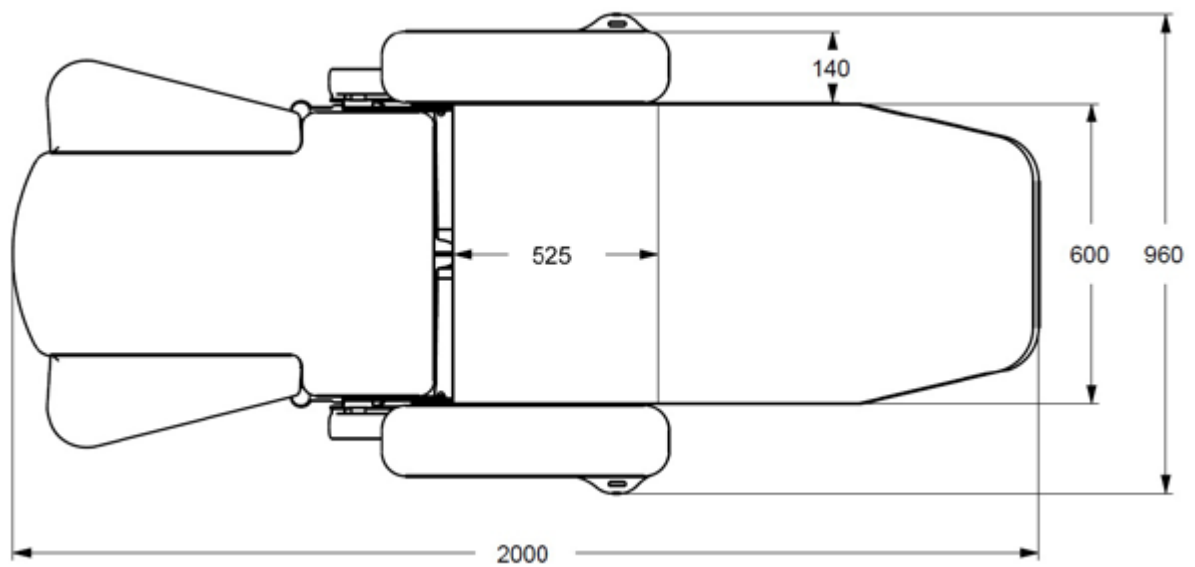



Data sheet

Codice CND / EMDN	V08020101
GMDN Code	38447
No. progressive/R	2433953/R
Basic UDI-DI	805771740POLPE
Product Code	GB0005.SP-2.0 Steel painted
Intended Use	Multifunctional armchair
Manufacturer	GARDHEN BILANCE SRL
Total weight	60 kg
Ruote	No. 4 stainless swiveling 100 mm with braking system
Power supply	100 - 240 Vac 50/60 Hz
Electric motors	Low voltage (24 V)
Max absorption	4A
Insulation	IP44
Maximum Scope	200 kg
Safe Workload	235 kg

Main Cord



Type F - Schuko



Type L - 16A



Type I - Australia



Type B - Canada and USA



Type G - UK plug



Gardhen Balance S.r.l.

Via G. Luraghi c/o Consorzio Il Sole - lot 5 - 80038 Pomigliano d'Arco (NA) - Italy

Tel. +39.081.8692160 pbx - Fax +39.081.8692460

E-mail: info@gardhenbalance.it

PEC: gardhenbalance@pec.it

www.gardhenbalance.it



Images reproduced in this sheet are for reference only and dimensional values are nominal. It is the manufacturer's right to make changes to the designs without prior notice. Refer to the accessories catalog and/or commercial proposals for the identification of standard equipment and the selection of any additional components. In case of discrepant data, the contents of the commercial proposals shall prevail.

Dep. EN_69 rev. 28/05/2024