Stephen H Comfort 6.0 Multifunctional chair



SEATING

Frequent Usage

Multifunctional Infusion therapy

Day Hospital

Blood Donation

Blood Sample

Emergency Room

Dialysis

Patient transport

Dermatology



HightCOMFORT Inpatient and outpatient therapy







Multi-adjustable armrests

Armrests ergonomic armrests equipped with special manual multiple adjustment mechanism for positioning in multiple directions. In addition, the support placed bilaterally on the armrests (equipment code AC25P) allows convenient storage of the hand control when not in use.



Sunchronized with the backrest



Swivable



Height adjustable



Removable



Tiltable









Movements

Independent back and leg sections


Independent movements of the back and leg sections are controlled by dedicated buttons on the hand control.





Movements Trendelenburg

The Trendelenburg position can be obtained at the dedicated button on the hand control.

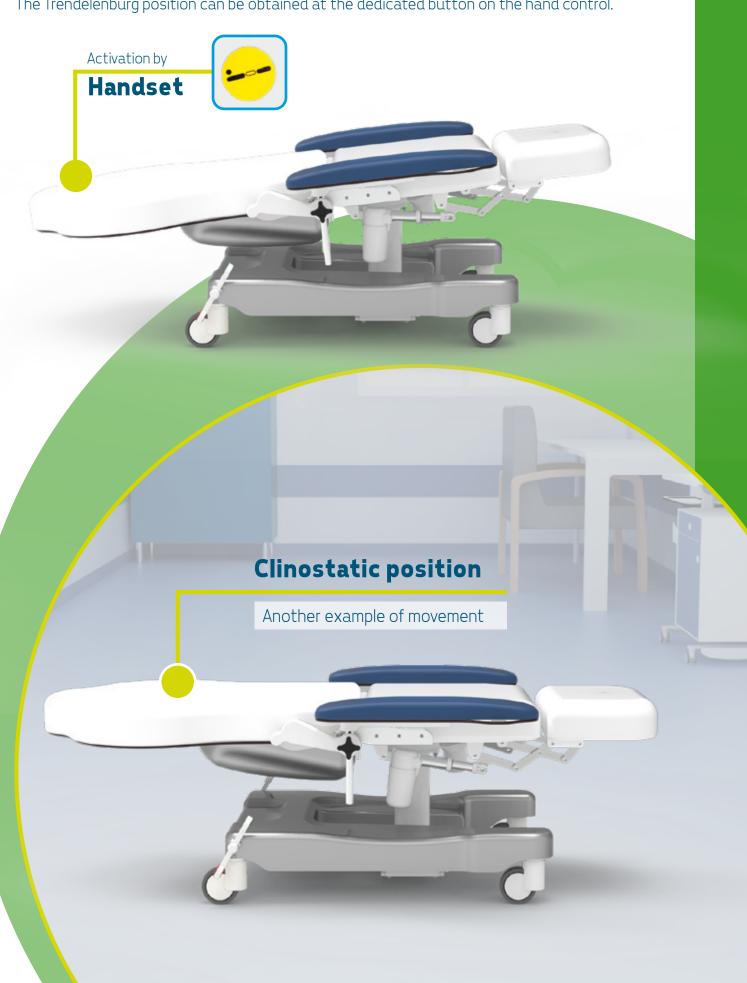



Other features of the medical device

- Antibacterial thermoset powder coated steel frame with wide assortment of colors.
- Top divided into 3 sections: back, seat and legs.
- By means of two dedicated motors, the back and leg sections can be adjusted independently and simultaneously Equipment code AC12.
- Leg section: 90° to the seat Equipment code AC63.
- The armrests are padded and upholstered in soft fabric. Thanks to a special clamp, they are: integral and synchronous to the movement of the backrest, swivel, height-adjustable, removable, tilt up/down Equipment code AC30/C1.
- Upholstery made of class 1IM leatherette, a material that can be sanitized with commonly used detergents and medical surgical aids, with antibacterial, antifungal, stain-resistant properties. Color customization upon request Equipment code AC87/1.
- No. 4 durable, antistatic, stainless, swivel wheels Ø 100 mm.
- Directional braking system useful for ensuring stability of the device-Equipment code AC7P1G.
- All chair movements are controlled by a perimeter hand control equipped with a 3.5 m coiled cable. The hand control is easy to use as it is made with large raised pictograms and is also suitable for visually impaired patients. The hand control is equipped with 10 keys with 2 memories. For patient safety, the hand control is equipped with a key for control inhibition Equipment code AC34P. IP66 insulation.
- Accommodation placed on the side on the armrests for push-button panel Equipment code AC25P.
- Cushions made of soft foam.
- ABS cover for the base Equipment code AC83
- ABS backrest cover Equipment code AC86C.

Do you want to change the color of the upholstery and frame? Frame or click the QR codes and find out how to customize your device by browsing our color cards*. *If no preference is expressed, the device will be supplied with the available color. *Upholstery Lethearette *Silicon *Frame Frame





Do you want to improve the device?

Frame the QR code and discover all the accessories available by browsing the catalog.

<u>F</u>

Technical Sheet

CND/EMDN Code **V08020101**GMDN Code **38447**

No. progressive/R 2430054/R

Basic UDI-DI **805771740POLPE**

ID Product GB0004.SP-6 Painted steel

Intended Use Multifunctional chair
Manufacturer GARDHEN BILANCE SRL

Total weight Standard 70 kg

Wheels No. 4 stainless swiveling 100 mm

with braking system

Power 100 - 240 Vac 50/60 Hz

Electric motors Low voltage (24 V)

Max absorption 4A Insulation IP44

Max load capacity Standard 200 kg Safe Workload Standard 235 kg

Main cord



Type F Schuko



Type L 16A



Type I Australia




Type B Canada and USA



Type G UK plug



Dimensions

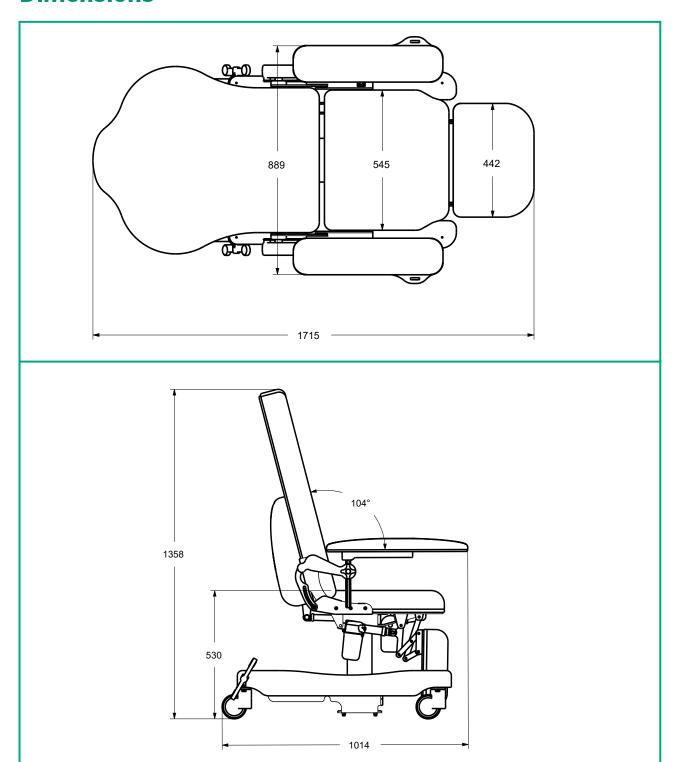

www.gardhenbilance.it











Via G. Luraghi c/o Consorzio II Sole - lotto S -80038 Pomigliano d'Arco (NA) - Italy Tel. +39.081.8692160 pbx - Fax +39.081.8692460 E-mail: info@gardhenbilance.it PEC: gardhenbilance@pec.it

